


Modular-Lab Software

General Overview

Modular-Lab Software

Intuitive and comprehensive software for programming and working with Modular-Lab synthesis equipment and quality control

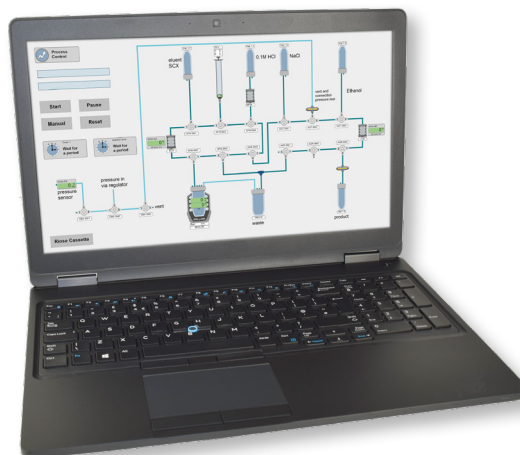

General Information

Modular-Lab Software is the renowned standard solution coming with all systems from the Modular-Lab family. Each synthesis device will be completely controlled and run by it. It comes in open as well as closed versions. Closed versions (as used with Modular-Lab easy) are non-editable and only support the use of pre-validated synthesis projects requiring a dedicated system setup. A wide range of according projects for different tracers is available via Eckert & Ziegler.

The open, editable software version (available for Modular-Lab Standard and PharmTracer) furthermore allows programming and creation of project templates for the development of synthesis processes. The software works with an intuitive graphical interface to create projects by drag and drop of symbols. Connected modules of the synthesis platform are automatically recognized and can be integrated into the process. An easy-to-use Macro Wizzard facilitates the creation of program flowcharts. The editable Modular-Lab Software versions require the use of a Dongle License in order to access and run the program.

Features

The Modular-Lab software supports users to comply with the highest standards of GMP, cGMP, GAMP 5 and 21 CFR part 11 regulations. An access control with four defined security levels meets all demands for the security of the process data and the application-specific setup. The logging of all system and user operations runs fully automatically.

Additional features include:

- Manual and automatic control of programmed processes
- Automated database download from PLC onto notebook after each project download
- Automatic backup function with adjustable settings to change intervals and limit storage usage
- System status messages
- Barcode reader
- Errors and other comments pertaining to specific process
- Display of steps in separate message window
- Four timers with countdown function
- Help function
- Trend viewer
- Calculation and comparison of mathematical and logical terms
- Service and maintenance information

Audit Trail and Reporting

The implemented audit trail function records all necessary information including process data (process steps, sensors), user data (user management, login, etc.), downloading of programs as well as occurring errors and messages.

The tracked data including time stamps is available as a comprehensive batch report (PDF file) after the end of a process.

Eckert & Ziegler Eurotop GmbH

Robert-Rössle-Strasse 10
13125 Berlin
Germany
Phone: +49 30 941084 197
Fax: +49 30 941084 470

medical@ezag.de
www.medical.ezag.com

Your contact in USA & Canada: Eckert & Ziegler Radiopharma, Inc.

25 Upton Drive
Wilmington, MA 01887
USA
Phone: +1 508 497 0060
Fax: +1 508 497 0061

InfoRAU@ezag.com
www.medical.ezag.com

